

LASAL

For a better experience in the water



September 2017

LASAL Flotation Cushions provide:

- Comfort & safety
- Stability in the water
- Maximum mobility
- Maximum movement
- Maximum joy

FLOAT



LASAL A/S, Industrivej 15, DK-8881 Thorsoe

E-mail: info@lasal.dk



+45 8696 6666



www.lasal.dk

VAT No. DK 1437 8677

Cushions for positioning in the pool can be used during training and different kinds of water therapy. They are developed specifically for disabled children and adults, but can also be used by people who want extra support in the water.

The cushions great float ability combined with skin friendly materials ensure the user nice and safe experience in the water.

With a good positioning in the pool the following is easily obtained:

- Increased mobility
- Increased stability
- Increased comfort and safety

For the user this means a more freely experience of being in the water.

For therapist or helper it improves the possibility of carrying out exercises or training that normally requires a second assistant.



LASAL Flotation Cushions are black on one side and blue on the other. We recommend to use the cushion with the black side against the skin as this side is the more smooth side.

Despite LASAL Flotation Cushions' great floatability they should always be distinguished from life jackets or other rescue equipment. Users of the Flotation Cushions must never stay in the water on their own or unattended and no longer than an arm's length away from assisting staff.

Extra attention should be paid to children below 3 years and people who are unable to swim or move on their own.

Collarfloat



9401, 9411, 9415, 9421

The Collarfloat is designed for supine or upright exercise in the water and prevents the user from getting his head under water.

Place the collar around the user's neck and close the front lengthwise of the chest. Next, tighten the strap around the cushion corresponding to just below the chin. This gives an upright position and the greatest possible freedom of movement in the water where only the head is above the water.

The long belt could be tightened around the back of the user in order to help the body up in a horizontal position.

The float cushion may advantageously be provided before the user enters into the water.



Banana-shaped Flotation Cushion



9101, 9111, 9811

Banana-shaped Flotation Cushions are available in more sizes and have very large floating effect. The cushions are typically placed in the armpit and around the user and is very useful for back, stomach and standing exercises.

For the more experienced swimmer this cushion may also be used for balance training and sitting exercises.



IMPORTANT

This cushion is not designed to be fastened around the body.
Learn more on page 9.

T-shaped Flotation Cushion



9017, 9027

The T-shaped Flotation Cushion is suitable for floating support of the lower limbs, either by placing the cushion behind the knees or the lower limbs.

For people with a certain stability in the water, the cushion could be used as floating support in a prone position. The cushion is placed in the armpits with T's foot pointing upwards towards the face. T's wings act as stabilizers in the water and T's foot as the main support. With this positioning the shoulders are freely mobile.




IMPORTANT

This cushion is not designed to be fastened around the body.

Learn more on page 9.

U-Cushion

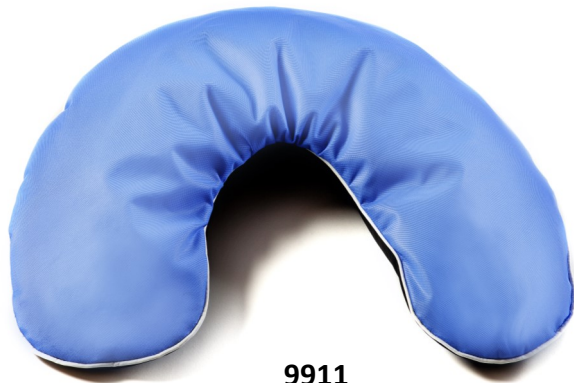

9911

The U-Cushion can be used by both adults and children and has several applications. Behind the neck in a supine position, under the abdominal or upper arms at breaststroke, or simply as extra flotation effect to the arms or the legs.



IMPORTANT

The U-Cushion is not designed to be fastened around the body. Learn more on page 9.

Neck Collar

The Neck Collar is fitted with grip tape closure, whereby the collar can be placed very close to the head and around the neck. The collar can be combined with other flotation cushions or used separately by those who wish or need just a little support during swim workouts.



9311, 9321

Bathing Flower



9600, 9610

The purpose of the Bathing Flower is similar to that of a bathing ring. The user can float leisurely in the flower, or be active in the supine as well as prone position. The petals and shape of the flower makes it easy for users to get their arms closer to the body. Depending on spasticity and movement pattern the flower can be used alone or combined with other LASAL Flotation Cushions. Available in two sizes.



Float Stabilizer



9620, 9630

The user is placed between the Float Stabilizer's two cushions in either sitting, standing or lying position. With the stabilizer the user can move actively around in the water in a supine or prone position, or use the cushion as extra support when floating.

For the more experienced user the stabilizer can be used for more active exercises or for balance training.

IMPORTANT

! The Stabilizer is not designed to be fastened around the body. Learn more on page 9. !

Specifications

Item No.	Product	Length	Width	Neckhole
9017	T-shaped Flotation Cushion, L	80 cm	54 cm	-
9027	T-shaped Flotation Cushion, S	63 cm	40 cm	-
9101	Banana-shaped Flotation Cushion, L	148 cm	38 cm	-
9111	Banana-shaped Flotation Cushion, M	122 cm	40 cm	-
9811	Banana-shaped Flotation Cushion, S	85 cm	40 cm	-
9311	Neck Collar, L	-	21 cm	24 cm
9321	Neck Collar, S	-	16 cm	24 cm
9401	Collarfloat, L	68 cm	47 cm	38 cm
9411	Collarfloat, M	65 cm	43 cm	32cm
9415	Collarfloat, S	59 cm	38,5 cm	26 cm
9421	Collarfloat, XS	54 cm	34 cm	26 cm
9911	U-Cushion	62 cm	21 cm	-
9600	Bathing Flower, S	-	74 cm	65 cm
9610	Bathing Flower, L	-	81 cm	89 cm
9620	Float Stabilizer, S	66 cm	19/24,5 cm	-
9630	Float Stabilizer, L	88 cm	24/29 cm	-

Description of Materials

Outer cover, facing up:	100% Polyester.
Outer cover, facing up:	100% Polyester.
Outer cover, facing down:	100% Polyester, mesh.
Filling:	Flame retarded EPS beads (Expanded Polystyrene).

Cleaning and maintenance

After use it is essential that the outer cover is rinsed thoroughly in non-chlorinated water. When drying the cushions must be placed with the blue mesh facing downwards in order to allow the water to drain from the cushions.

The cushions take quite a long time to become completely dry and should therefore be used exclusively in "wet" environments.

Disposal

The cushions can be disposed of together with common household waste.



IMPORTANT INFORMATION

These cushions are not designed to be fastened around the body and it is therefore necessary to ensure that they do not float away from the user. Either user or helper must be able to hold on to the cushion in order to avoid getting the head under water.

LASAL Flotation Cushions are exclusively intended for use in chlorinated water.



!!! ATTENTION !!!

Avoid contact with sharp and pointed items/edges and rough surfaces against the cushion



We recommend always to read the manual before using the product.

The responsibility for the safety of the bathers, the risk of drowning or other injury always rest with the assistant staff whether LASAL Floating Cushions have been used or not.



LASAL A/S, Industrivej 15, DK-8881 Thorsoe

E-mail: info@lasal.dk



+45 8696 6666



www.lasal.dk

VAT No. DK 1437 8677